

MATTRESSES FOR CENTRELLA® BEDS

Our mattresses and covers are a perfect fit for your Centrella beds. The standard model has two foam layers, while the Reactive Air mattress adds extra comfort and protection against skin breakdown. Both feature top-layer memory foam, and the Reactive Air includes air cells that respond to movement — no electricity or loud motors needed.



WATERFALL ZIPPER FLAP
The waterfall zipper flap covers the zipper, so fluids pass over it without absorbing into the mattress.

PERFECT FIT
All our mattresses are accurately made to be a perfect fit on your beds.

POLY-C
Our mattresses feature a durable polyester coated knit with a polycarbonate additive. This material provides greater resistance to aggressive cleaners, including bleach and hydrogen peroxide.

EXPANDABLE HEEL SECTION
Provides consistent support and pressure redistribution to aid in the prevention of pressure injuries.

GEL-INFUSED MEMORY FOAM
Memory foam reduces pressure to body regions that are susceptible to pressure injuries. Gel particles embedded in the foam help to draw heat away from the body while offering greater structural integrity to the foam.

AIR CELLS
Our Reactive Air mattresses contain a series of interconnected air cells that deliver dynamic response to patient movement without electricity or loud motors.

5 YEAR
FIVE YEAR WARRANTY
We stand behind all our bed mattresses with a five year warranty.

FULL-PERIMETER ZIPPER
Our full-perimeter zippers provide for complete separation of top and bottom covers. This makes it easy to replace a worn top cover while keeping the bottom cover intact.

RF WELDED SEAMS
Rather than a sewn seam, our seams are sealed using radio frequency welding to prevent fluid ingress.


MULTI-DIRECTIONAL STRETCH
Our breathable four-way stretch fabric provides added comfort and durability. It also reduces friction for easier patient repositioning.

FIRE SAFETY
All of our mattresses meet Fed 1633 fire safety standards. Our full-coverage flame barriers, standard on all mattress assemblies, provide additional protection and durability.


THOROUGH TESTING
All of our mattresses are designed, manufactured, and tested to the S3I standard. For more information about the Support Surface Standards Initiative, please visit <https://npiap.com/page/S3I>.

MADE IN THE USA
All of our mattresses are manufactured in the USA.

FOAM MATTRESS CONSTRUCTION

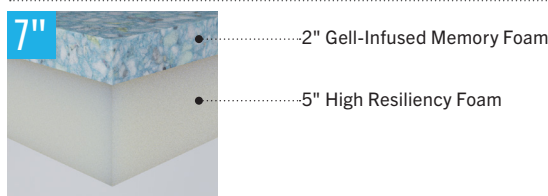

REACTIVE AIR MATTRESS CONSTRUCTION

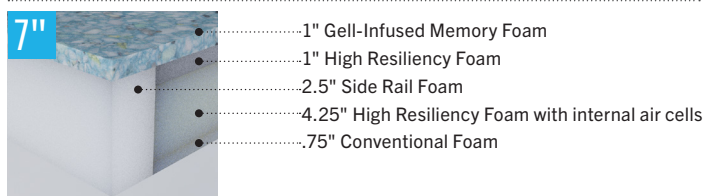

Type	Width	Length	Thickness	Patient Weight Limit	FULL MATTRESS	TOP COVER	TOP AND BOTTOM COVER
					Part Number	Part Number	Part Number
Foam	36"	88"	7"	550 pounds	29707	29709	29743
Reactive Air	36"	88"	7"	550 pounds	29708	29709	29743
Max Surface Top Cover	36"	88"	—	—	—	29680	—

Centrella is a registered trademark of Hillrom

©2025 Curbell Medical Products, Inc. MAP1922A